



SEC Separation of Antibiotics Biapenem on Zenix™-100(7830)

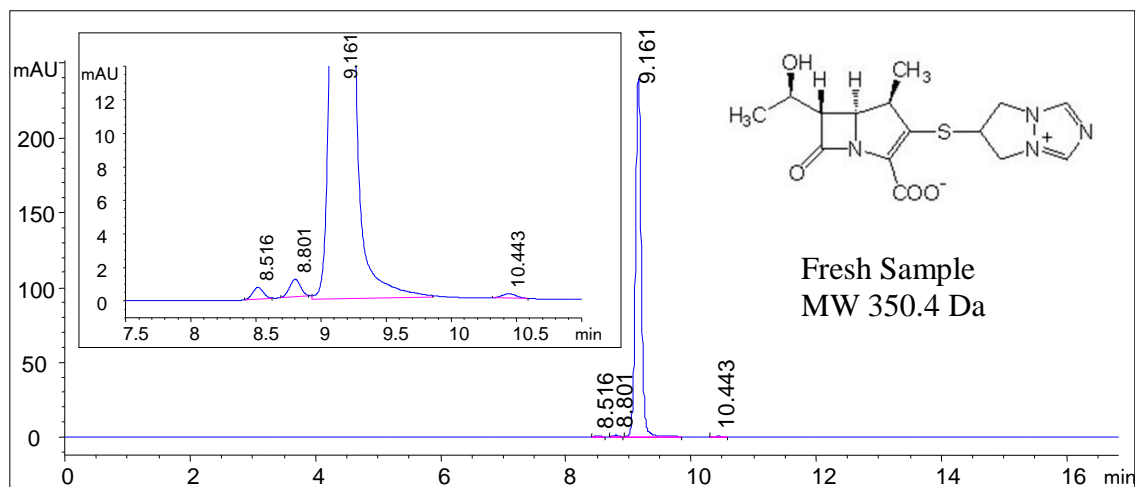

Columns:	Zenix™-100 (3 μm, 100 Å, 7.8x300 mm)
Mobile Phase:	150 mM sodium phosphate buffer, pH 7.0
Flow Rate:	1.0 mL/min
Temperature:	Ambient
Detection:	UV 215 nm
Injection volume:	20 μL
Samples:	Biapenem (0.2 mg/mL in mobile phase)

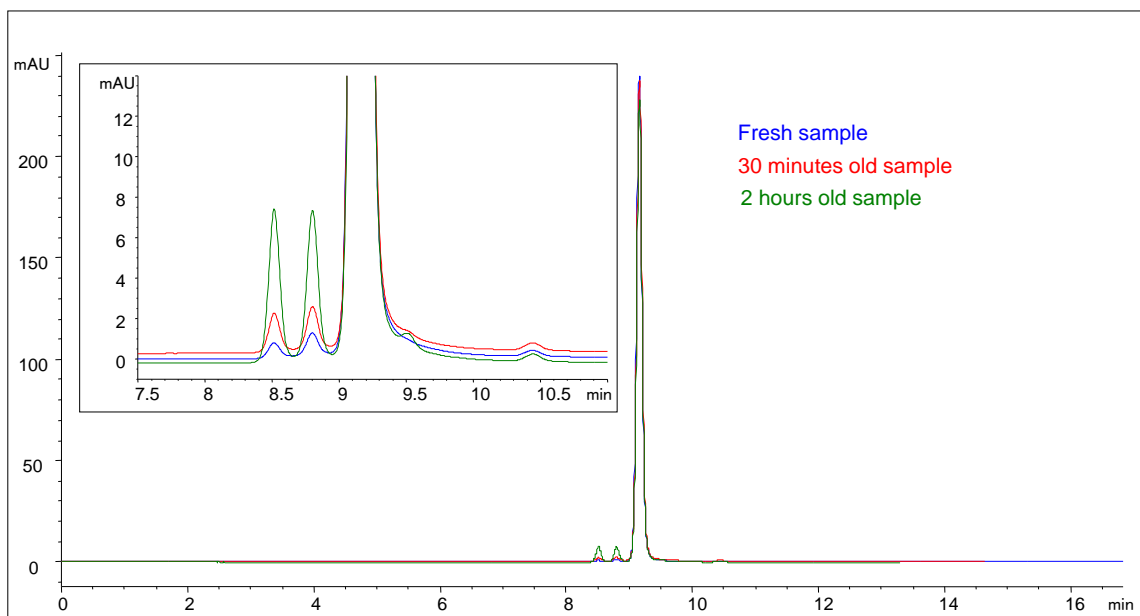

Keywords: Size exclusion, Zenix, Antibiotic, high resolution, Biapenem, Pharmaceuticals