




Zenix™ Column Lot-to-Lot Reproducibility for MAb

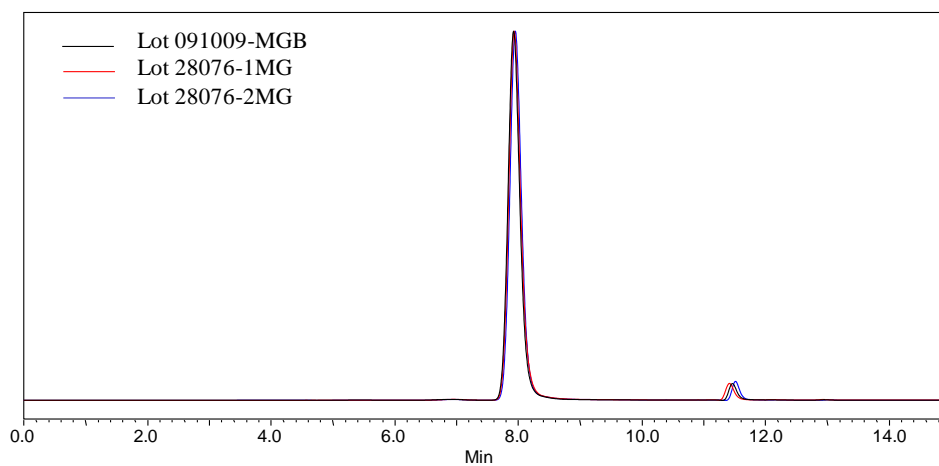

Column: Zenix™ SEC-300 (3 μ m, 300 Å, 7.8x300 mm)
Mobile Phase: 150 mM Phosphate Buffer, pH 7
Flow Rate: 1.0 mL/min
Temperature: Ambient
Detection: UV 214nm
Concentration: 2.5 mg/ml
Injection Volume: 10 μ L
Sample: MAb

Keywords: Size exclusion, Zenix, monoclonal antibody, lot-to-lot